MEDIPHOT 937/943 X-RAY

COLENTA Medical X-Ray Film Processors: MP937 and MP943


For all standard X-Ray Films - Certified for Mammography - Film processing.

MP937/943: max. 37 (43) cm **DEV-Time** 18sec to 99 sec

Cycle Time 1'24"min to 7'43"min

Intake Speed 106cm/min to (min.) 20cm/min

(max.)

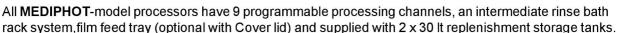



Standard Features and accessories:

- compact / robust design
- open stand (optional covered)
- 2 rep. storage tanks 30L (for Dev and Fix)
- optical sensorbar for touchless measuring incoming film area

VERSION I Media Exit on **Top Cover Lid** **VERSION II Rear Media Exit** in Receiving Basket

Certified

Company








The best results in "Low Volume" Tank designs

All Colenta Processors are equipped with precision temperature control of all processing chemistry with automatic Developer Tank cooling to prevent overheating. Replenishment is preheated by injection into the inline chemical heating systems, which are controlled to pre-set programmed temperatures within +/- 0,2 °C.

TRANSPORT SYSTEM:

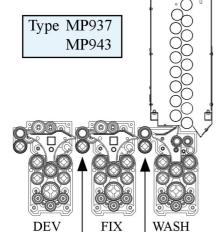

Rinsed cross over roller







Microprocessor controlled Displa

Optional accessories:

Rack Carrier Tray/Drip Tray for easy handling the transport racks.

Cover Lid for Film Feed Tray





Closed stand for processor



FUTUREMED İÇ VE DIŞ TİCARET LTD. ŞTI

Tel: +90 212 210 2142

GSM: +90 533 784 4295

Web site: www.futuremed-tr.com/tr

Email: info@futuremed-tr.com

document no.: X1-5 1.pdf